



Technology for Life Science

Kowa Applanation Tonometer



HA-2



Convenient Item For Daily Clinical Practice



● Standard set of HA-2

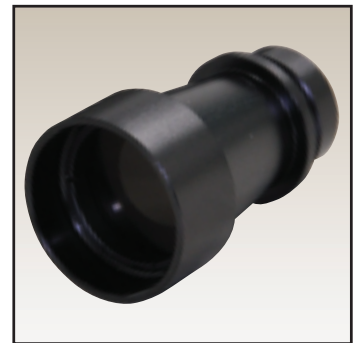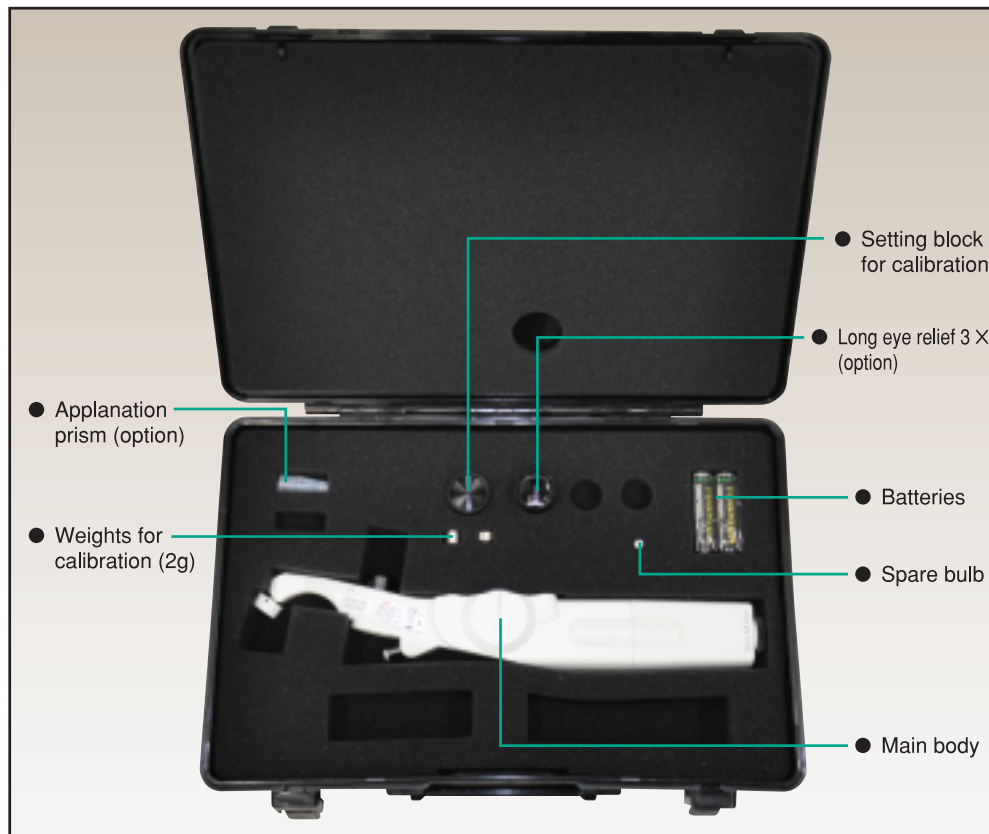

Long eye relief 3 X (option)

● Main Features

- ◆ Compact, lightweight and well-balanced; easy to handle with one hand.
- ◆ Unique one-spring mechanism assures correct reading regardless of patient posture.
- ◆ Built-in illumination permits you to use the instrument regardless of the lighting conditions of consultation room.
- ◆ Direct reading type precision dial requires no conversion table to change the readings.
- ◆ Accessories such as adjustable head rest and interchangeable eyepiece provide most suitable conditions for measurement.

● Specifications

Type	Hand-held applanation tonometer
Applanation prism	Doubling prism interchangeable (option)
Scale range	0~60 mm Hg (1 mm Hg division)
Illumination	
a) ON/OFF	Interlocked with thumb wheel
b) Bulb	Lens bulb 2.5 V, 0.25 A interchangeable
c) Battery used	1.5 V, 2 pcs Interchangeable
Eyepiece	
Magnification	6 X
Head rest	Movable
Length	290 mm
Weight	240g/0.5lbs

* All specifications and external design are subject to modification of change without prior notice.



Distribution name : Kowa Applanation Tonometer HA-2



Kowa Company, Ltd.

World Sales Headquarters

4-14, Nihonbashi-honcho 3-chome, Chuo-ku, Tokyo 103-8433 Japan
 Phone: 81(3)3279-7639 Facsimile: 81(3)3279-7541
 URL: <http://www.kowa.co.jp/e-life/>

Hamamatsu Factory

3-1, Shinmiyakoda 1-chome, Kita-ku, Hamamatsu City, Shizuoka Pref., 431-2103 Japan



Kowa Europe GmbH

Immermannstrasse 43B
 40210 Duesseldorf, F.R. Germany
 Phone: 49(211)1793540
 Facsimile: 49(211)161952
 URL: <http://kowa-europe.com/>

Kowa Optimed, Inc.

20001 S. Vermont Ave.
 Torrance, CA 90502, U.S.A.
 Phone: 1(310)327-1913
 Facsimile: 1(310)327-4177
 URL: <http://kowa-usa.com/>



ATO6_1097200MA/K
 Printed in Japan.